



www.medi-plinth.co.uk



PATIENT HANDLING BROCHURE

Designed for comfort and function

Stalham Patient Handling

CONTENTS

STALHAM 150 Lifting Hoist	3
STALHAM 150 Sit-stand Lift	4
STALHAM 20 Portable Commode	7
STALHAM 64 Portable Bariatric Commode	9
STALHAM TILT Portable Shower Commode	10
STALHAM COMPACT Shower Transfer Chair	12
STALHAM BRAKING Shower Transfer Chair	13
STALHAM STABLE Commode/Shower Chairs	15

Patient Hoists

STALHAM 150 LIFTING HOIST



Structure

Epoxy coated steel frame.

Complete lifting hoist can be assembled in just a few minutes.

Utilises electric lifting arm.



Comfort

A range of comfortable slings are available to fit patient requirements.

Ergonomic handles.



Ease Of Handling

Wired remote control for user controls in any position around the lift.

Capable of lowering patient to the ground.

The base frame has a very low construction to allow for very good under-ride possibilities.

Undercarriage is expandable as needed by means of foot pedal.

Mounted on four swivelling castors for portability. Each rear castor is braked.



Safety

Features easy to access emergency stop button.

Lift can allow for emergency mechanical lowering if required.

Additional battery backup.

Medical grade electric lifting arm.

Patient Hoists

STALHAM 150 SIT-STAND LIFT



Structure

Epoxy coated steel frame.

Complete lifting hoist can be assembled in just a few minutes.

Utilises electric lifting arm.



Comfort

Leg straps ensure comfort and stability for patient when lifting.

180° grab handles support patient.

Angled foot rests secure feet in place when in use.



Ease Of Handling

Wired remote control for user controls in any position around the lift.

The base frame has a very low construction to allow for very good under-ride possibilities.

Undercarriage is expandable as needed by means of foot pedal.

Mounted on four swivelling castors for portability. Each rear castor is braked.



Safety

Features easy to access emergency stop button.

Lift can allow for emergency mechanical lowering if required.

Additional battery backup.

Medical grade electric lifting arm.

Specifications

STALHAM 150 LIFTING HOIST

Hoist: 150 LIFTING HOIST

Lift height (min.)	55 cm
Lift height (max.)	174 cm
Underpin height	14 cm
Device weight	43 kg
Max. user weight	150 kg

Accessories:

Toilet sling (optional velcro fastener version)
Approved up to 250 kg weight



Universal sling (S, M/L, XL/XXL)
Approved up to 250 kg weight



Universal sling with headrest (M/L)
Approved up to 300 kg weight



Bariatric sling
Approved up to 300 kg weight




Net/bath sling (S, M/L, XL/XXL)
Approved up to 250 kg weight




Wide Hanger
4-point spreader bar 56 x 33 x 14cm



Horizontal sling (optional split leg)
Approved up to 250 kg weight

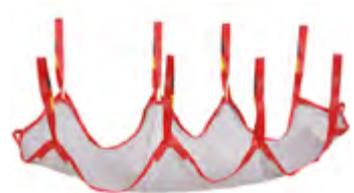

Hoist Scale PW500
Approved to class III, MPG class I

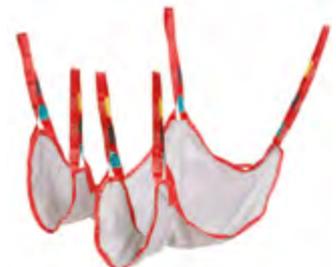

Specifications

STALHAM 150 SIT-STAND HOIST

Hoist: 150 SIT-STAND LIFT

Lift height (min.)	84 cm
Lift height (max.)	150 cm
Underpin height	14 cm
Device weight	53 kg
Max. user weight	150 kg

Accessories:

Back support sling

Approved up to 250 kg weight



Calf support belt



www.medi-plinth.co.uk

7-11 Holywells Road, Ipswich, Suffolk IP3 0DL – Call: +44 (0) 1473 212010

Portable Commodes

STALHAM 20



Structure

Utilises a completely synthetic frame making it durable and lightweight. Coated tube frame design allows for a very high resistance to wear

Strong, rigid materials also have the added benefit of shock dissipation.

Includes a moulded removable potty as standard (stainless steel rack also available)

Comfort

Moulded backrest is ergonomically designed to be comfortable for users.

Available as fixed, adjustable or wide models. Offering adjustable height or increased width to accommodate patients of different sizes.

Footrests are height adjustable and swing away when not required.

Moulded synthetic rubber over-seat which securely encloses the aperture and is kind to patient skin.

Ease Of Handling

Push handle integrated in back frame for easy and accurate portering.

Fitted with four Ø 125mm braked castors

Arm rests which are easily retractable to enable secure transfer and patient access.

Rubber over-seat has no seams and will not rip/ tear.

Easy Cleaning

Chair can be dismantled within 20 seconds without tools to allow for thorough cleaning.

Designed without welding to avoid rust or micro-fracture weak points. Seamless joints ensure dirt or bodily fluids cannot be trapped.

Capable of withstanding the most rigorous cleaning regime. Suitable for cleaning with NHS approved cleaning solutions and Clinell wipes (or other alcohol based wipes).

Specifications

STALHAM 20 RANGE

Commode: STALHAM 20

Overall dimensions	102 x 57 x 60 cm
Aperture seat height	51 cm
Under aperture height	49 cm
Over-seat height	53 cm
Seat width	47 cm
Seat depth	40 cm
Armrest height	26 cm
Backrest height	52 cm
Clears toilet of height	49 cm
Safe working load	130 kg
Product weight	15 kg

Commode: ADJUSTABLE HEIGHT MODEL

Identical dimensions but with height adjustment:

Height	101–107 cm
Aperture seat height	48–53 cm
Clears toilet of height	47–520 cm

Commode: WIDE MODEL

A wider design with larger safe working load:

Overall dimensions	102 x 62 x 60 cm
Seat width	52 cm
Seat depth	45 cm
Safe working load	150 kg
Product weight	16 kg

Portable Commodes

STALHAM 64 BARIATRIC



Structure

Utilises a completely synthetic frame making it durable and lightweight. Coated tube frame design allows for a very high resistance to wear

Strong, rigid materials also have the added benefit of shock dissipation.

Heavy-duty backrest with aluminium construction to secure patient.

Includes a bucket device and bucket.

Comfort

Moulded back section features contoured rubber backrest for optimum patient comfort.

Wider seat is designed to accommodate larger patients.

The foot rests are height adjustable and swing away when not required.

Moulded synthetic rubber over-seat which securely encloses the aperture and is kind to patient skin.

Ease Of Handling

Arched rear frame section acts as ergonomic push handle and allows for easy and accurate portering.

Fitted with four Ø 125mm braked castors

Arm rests which are easily retractable to enable secure transfer and patient access.

Rubber over-seat has no seams and will not rip/ tear.

Easy Cleaning

Chair can be dismantled within 20 seconds without tools to allow for thorough cleaning.

Designed without welding to avoid rust or micro-fracture weak points. Seamless joints ensure dirt or bodily fluids cannot be trapped. Seamless joints ensure dirt or bodily fluids cannot be trapped.

Capable of withstanding the most rigorous cleaning regime. Suitable for cleaning with NHS approved cleaning solutions and Clinell wipes (or other alcohol based wipes).

Portable Commodes

STALHAM TILT SHOWER COMMODE



Structure

Utilises a completely synthetic frame making it durable and lightweight. Coated tube frame design allows for a very high resistance to wear

Strong, rigid materials also have the added benefit of shock dissipation.

Available with electric or hydraulic lift mechanism.

Chair frame extends to recline seat and backrest.

Also includes a bucket device and bucket.



Comfort

Comfortable, adjustable head rest which combined with the tilting function allows for endless positions with this versatile chair.

Sturdy, contoured back, head and leg supports have been ergonomically designed for optimum patient comfort.

The foot rests are height adjustable and swing away when not required.

Moulded synthetic rubber over-seat which sits over the aperture securely and is kind to the skin.



Ease Of Handling

Easy to use level mechanism allows fine adjustment of the chair back rest angle.

A push handle is integrated for easy and accurate portering.

Headrest tilt and height are adjusted with simple to use clip mechanisms.

Fitted with four Ø 125mm braked castors

Arm rests which are easily retractable to enable secure transfer and patient access.

Rubber over-seat has no seams and will not rip/tear.



Easy Cleaning

Chair can be dismantled within 20 seconds without tools to allow for thorough cleaning.

Designed without welding to avoid rust or micro-fracture weak points. Seamless joints ensure dirt or bodily fluids cannot be trapped. Seamless joints ensure dirt or bodily fluids cannot be trapped.

Capable of withstanding the most rigorous cleaning regime. Suitable for cleaning with NHS approved cleaning solutions and Clinell wipes (or other alcohol based wipes).

Specifications

STALHAM 64 BARIATRIC & TILT SHOWER COMMODORE

Commode: 64 BARIATRIC

Overall dimensions	104 x 68 x 90 cm
Seat height	56 cm
Seat depth	48 cm
Seat width	64 cm
Toilet height with bucket	44 cm
Toilet height w/o bucket	47 cm
Safe working load	200 kg
Product weight	21 kg



Dialysis Chair: TILT SHOWER COMMODORE

Overall dimensions	120 x 111 x 65 cm
Seat width	54 cm
Seat depth	43 cm
Length when extended	144cm
Armrest height	76 cm
Push handle height	104 cm
Backrest tilt from seat	103°–128°
Headrest Adjustment	± 12 cm
Headrest Rotation	360°
Safe working load	150 kg
Product weight	29 kg

Shower Transfer Chairs

STALHAM COMPACT



Structure

Utilises a completely synthetic frame making it durable and lightweight. Coated tube frame design allows for a very high resistance to wear

Strong, rigid materials also have the added benefit of shock dissipation.

Compact design for space-saving benefits.

Includes a moulded removable potty as standard (stainless steel rack also available)

Comfort

Moulded backrest is ergonomically designed to be comfortable for users.

Footrests are height adjustable and swing away when not required.

Seat pad is contoured for patient comfort and features a rear drainage channel for water run-off

Ease Of Handling

Push handle integrated in back frame for easy and accurate portering.

Fitted with four Ø 125mm braked castors

Arm rests which are easily retractable to enable secure transfer and patient access.

Rubber over-seat has no seams and will not rip/ tear.

Easy Cleaning

Chair can be dismantled within 20 seconds without tools to allow for thorough cleaning.

Designed without welding to avoid rust or micro-fracture weak points. Seamless joints ensure dirt or bodily fluids cannot be trapped.

Capable of withstanding the most rigorous cleaning regime. Suitable for cleaning with NHS approved cleaning solutions and Clinell wipes (or other alcohol based wipes).

Shower Transfer Chairs

STALHAM BRAKING



Structure

Utilises a completely synthetic frame making it durable and lightweight. Coated tube frame design allows for a very high resistance to wear

Strong, rigid materials also have the added benefit of shock dissipation.

Heavy-duty backrest with aluminium construction to secure patient.

Accessory points are available for the addition of: IV pole, safety belt and arm pads



Comfort

Moulded back section features contoured rubber backrest for optimum patient comfort.

The foot rests are height adjustable and swing away when not required.

Seat pad is contoured for patient comfort and features a rear drainage channel for water run-off



Ease Of Handling

Arched rear frame section acts as ergonomic push handle and allows for easy and accurate portering.

Fitted with four Ø 125mm castors, two rear locking which can either swivel freely, lock in-line or fully lock.

Arm rests which are easily retractable to enable secure transfer and patient access.

Rubber over-seat has no seams and will not rip/ tear.



Easy Cleaning

Chair can be dismantled within 20 seconds without tools to allow for thorough cleaning.

Designed without welding to avoid rust or micro-fracture weak points. Seamless joints ensure dirt or bodily fluids cannot be trapped.

Capable of withstanding the most rigorous cleaning regime. Suitable for cleaning with NHS approved cleaning solutions and Clinell wipes (or other alcohol based wipes).

Specifications

STALHAM COMPACT & BRAKING

Shower Chair: COMPACT

Overall dimensions	102 x 57 x 60 cm
Seat height	53 cm
Seat depth	40 cm
Seat width	47 cm
Armrest height	26 cm
Backrest height	52 cm
Safe working load	130 kg
Product weight	15 kg

Shower Chair: BRAKING

Overall dimensions	102 x 64 x 76 cm
Seat height	56 cm
Seat depth	46 cm
Seat width	54 cm
Armrest height	26 cm
Backrest height	52 cm
Safe working load	130 kg
Product weight	15 kg

Shower Chairs

STALHAM STABLE COMMODE CHAIR & TOILET FRAME



Features

Manufactured with plastic and aluminium.

Designed to fit over the toilet.

Adjustable height with seven positions.

Approved up to 130 kg body weight.

Available with or without backrest.

Accessories

- Bucket
- Seat upholstery
- Travel carry bag
- Holder for disposable bags
- Disposable bags

Chair: STABLE COMMODE CHAIR

Seat height	42 – 57 cm
Seat width	46 cm
Seat depth	40 cm
Max required space	51 cm x 59 cm
Weight w/o accessories	~5 kg
Max. userweight	130 kg

Frame: STABLE TOILET FRAME & SEAT

Max required space	51 cm x 59 cm
Weight w/o accessories	~5 kg

Shower Chairs

STALHAM STABLE SHOWER CHAIR & STOOL



Features

Manufactured with plastic and aluminium.

Designed for patient showering.

Adjustable height with seven positions.

Approved up to 130 kg body weight.

Available with or without backrest.

Accessories

- Travel carry bag

Chair: STABLE SHOWER CHAIR

Seat height	42 – 57 cm
Seat width	40 cm
Seat depth	40 cm
Max required space	51 cm x 53 cm
Weight w/o accessories	~5 kg
Max. userweight	130 kg

Stool: STABLE SHOWER STOOL

Max required space	45 cm x 53 cm
Weight w/o accessories	~4 kg



www.medi-plinth.co.uk

7-11 Holywells Road, Ipswich, Suffolk IP3 0DL – Call: +44 (0) 1473 212010