
Medi Plinth

www.medi-plinth.co.uk



STALHAM BATHE 17

Height adjustable paediatric bath

Robust, space-saving option for children

Adjustable Bath

STALHAM BATHE 17



Structure

The bath tub is manufactured from fibreglass reinforced polyester. The frame and skirts are made of powder coated steel.

Attractive design, colour options for feet and skirts.

The smaller, space-saving design is suitable for environments with limited space, without sacrificing function, efficiency or comfort.

Lift system supports weights of up to 150 kg.

Comfort

Padded mattress included for patient comfort, is easily removable for cleaning.

Height adjustable from 66-96cm

Dual-synchronised lifting columns allow for quick and quiet height adjustment.

Factory-fitted thermostatic mixer regulates both the bathtub filling and the shower water temperature, to ensure a comfortable temperature for the user.

Ease Of Handling

Overflow control and pop-up drainage with hose. 40mm drainage can be attached to the floor or wall.

A pop-up drainage closes and opens by turning the overflow valve. No hands need to be in the water.

Adjustable feet increase safety and comfort if the floor should be uneven.

Shower handle has a soft rubber outlet nozzle which is easy to clean and that gives comfortable soft water dispersion.

Clearance below the tub allows the use of many different types of hoists and easy floor cleaning

Controls

Height adjustment is operated via a simple hand controller, with convenient placement hole.

The hand control can be operated from either side of the bathtub as required.

Specifications



STALHAM BATHE 17

Paediatric Bath: BATHE 17

Length	168.5 cm
Width	69.5 cm
Height	66–96 cm
Lift speed (no load)	16 mm/s
Safe working load	150 kg
Weight (empty)	230kg

Electric features:

Power Supply	230 V, 50 Hz
--------------	--------------

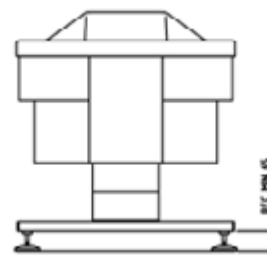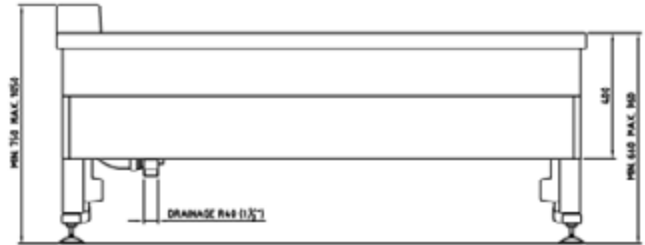

www.medi-plinth.co.uk

7–11 Holywells Road, Ipswich, Suffolk IP3 0DL – Call: +44 (0) 1473 212010